


THE APPLICATION


- Intuitive and easy to use
- Soon available on Google Play*
- Allows to monitor the patient's drainage data remotely



- Data export for further analyses

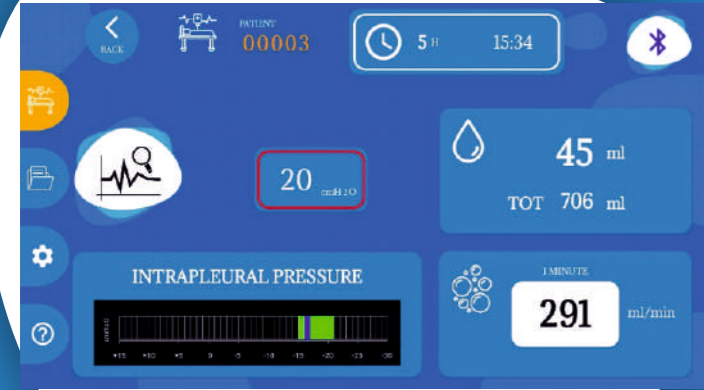

- Visualization of real-time graphs for every measurement available.

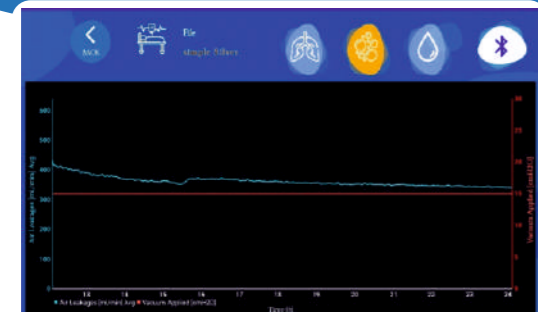

ACCESSORIES



code 10539
Chest unit protection filter
Packaging: 36 units for box



code 10547
Cap for Chest Drain
Packaging: 50 units for box



SilverLine app
Download and data processing software
COMING SOON


DRENTECH™ SILVERLINE



The TRULY EVOLVED drainage system

TECHNICAL FEATURES

- OLED display to enhance visibility in every situation
- Integrated anti-tilting sensor

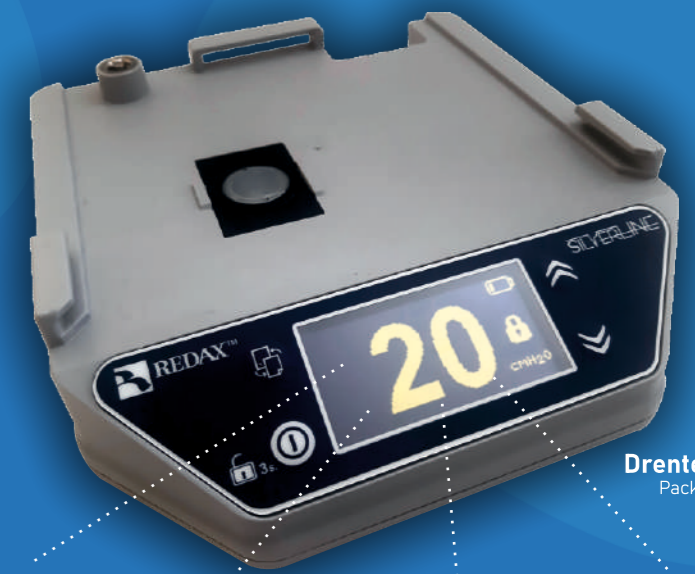

DrenTech SilverLine
Packaging: 1 unit for box
code 10198

Negative pressure adjustable from 0 up to -60cmH₂O




Air leaks monitoring



Liquid leaks monitoring



Trend of intrapleural pressure

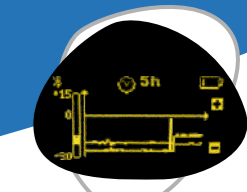

- Bluetooth connectivity to allow a remote monitoring of the patient
- Up to 150 hours of registration and data storage
- Available with 3 different disposable



REDAX™
MEDICAL DEVICES
Via Galileo Galilei, 18 - 46025 Poggio Rusco (MN) ITALY
Tel. +39 0386 830582 | Fax +39 0386 51898
E-mail: info@redax.it | www.redax.it

REDAX™
MEDICAL DEVICES

* "Google Play and the Google Play logo are trademarks of Google LLC."

DRENTECH™ SIMPLY

Hybrid single bottle drainage system

- 1200 ml collection chamber with integrated underwater seal valve
- Fast lock-in connector for patient tube



All-time clear visibility for air leakage

DOUBLE LUMEN


Allows a real measurement of intrapleural pressure



Double Lumen Tube_ single connector
Packaging: 24 units for box code **10541**



Double Lumen Tube_ Y connector
Packaging: 24 units for box code **10540**

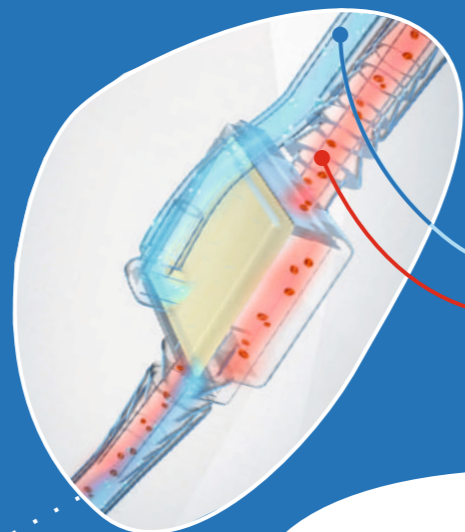

Hydrophobic membrane for the separation of AIR from LIQUID

- No need of milking or manipulation of the tube
- Less workload for hospital staff

BENEFITS

↑ **PATENCY** of tube > Faster lungs re-expansion

↓ **EFFORT** of the patient > Early mobilization of the patient

↓ **RESISTANCE** caused by liquid >

DRENTECH™ min SIMPLY

Hybrid single bottle drainage system

- 400 ml collection chamber with integrated underwater seal valve
- Fast lock-in connector for patient tube
- Designed for paediatric use



Dretech miniSimply_ bottle only code **10176**
Packaging: 12 units for box

DRENTECH™ SIMPLE SILVER

Two chambers drainage system

- 2000 ml collection chamber separated from underwater seal valve
- Fast lock-in connector for patient tube
- Designed for SilverLine use
- Collection chamber: a unique chamber to have a precise reading of liquid leakage



Dretech SimpleSilver_ bottle only code **10121**
Packaging: 4 units for box

Dretech Simply_ bottle only code **10171**
Packaging: 12 units for box